




VET LIFT STAINLESS

Veterinary Surgery and Treatment Table

- Made of high-quality stainless steel
- Electrochemically polished surface
 - > High corrosion resistance
 - > Plain and polished surface
 - > Optimal cleanability
 - > Strongly reduced affinity to coating
 - > Metallic pureness and chemical passiveness
- Electrical (230V power cord) or battery-powered height adjustment by foot switch from 31-101 cm
- Frame of maximum stability
- Plate (stainless steel or synthetic top) with edge and drain (130 x 60 cm) as well as 14 tie-down hooks
- Optionally fixed with 2 castors installed on one side (allow smooth movement of the table by lifting on the opposite side) or moveable by 4 double castors with blocking mechanism
- Top optionally fixed or tiltable (continuously adjustable tilt angle of the top (up to 11°) with comfortable one-hand operation supported by gas spring)

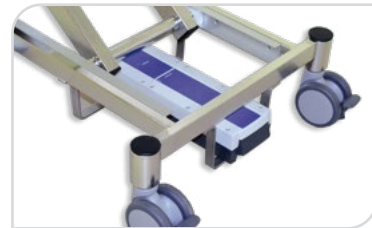

VET LIFT Stainless Steel

fixed with small castors on one side



Electrical height adjustable (power cord)

Top	Item Code
fixed stainless steel top	659.605
tilting stainless steel top	659.606
fixed synthetic top	659.615
tilting synthetic top	659.616

Battery-powered height adjustable

Top	Item Code
fixed stainless steel top	679.605
tilting stainless steel top	679.606
fixed synthetic top	679.615
tilting synthetic top	679.616

VET LIFT Stainless Steel

movable by 4 double castors (with blocking mechanism)



Electrical height adjustable (power cord)

Top	Item Code
fixed stainless steel top	659.601
tilting stainless steel top	659.602
fixed synthetic top	659.611
tilting synthetic top	659.612

Battery-powered height adjustable

Top	Item Code
fixed stainless steel top	679.601
tilting stainless steel top	679.602
fixed synthetic top	679.611
tilting synthetic top	679.612

